

AIRSON[®]

Air scaler



TeKne Dental


AIRSON®

For fast and simple tartar removal

Ease of installation

The instrument can be used on any unit, in place of the turbine handpiece, without requiring any additional adjustment.

It is available with the 3-hole, 4-hole or quick connection which, together with the MULTIWAY quick couplings, allow the instrument to rotate 360 degrees.

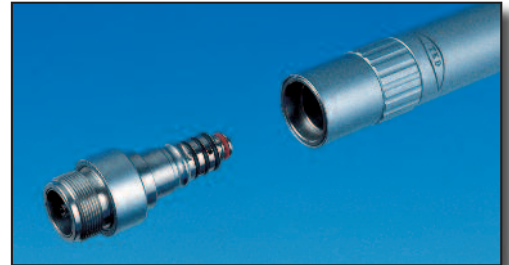

Power adjustment

The power of the vibrations can be adjusted with a rotating ring, making the prophylaxis more flexible and accurate.

Maintenance is extremely easy: a few drops of lubricating oil are sufficient for many days of use.




Inserts

Inserts available are:

- #1 Universal
- #2 Extra-fine
- #4 Spatula

Other patterns are available individually:

- #5 Amalgam condensation (small head)
- #6 Amalgam condensation (big head)
- #7 Crown and pivot loosening

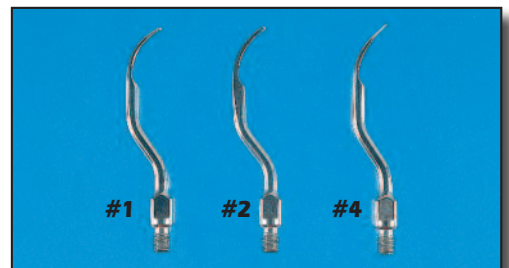

Insert replacement

The fitting of inserts is easy and secure with the supplied autoclavable key, which also protects against cross infection and injuries.



Technical data

Reference number:	REF 232.x0
Operating frequency:	6,000 ÷ 6,700 Hz
Recommended air pressure:	2.0 ÷ 2.8 bar (29 ÷ 41 psi)
Air consumption:	18 NI/min (@ 2.8 bar)
Weight:	65 g
Sterilization:	By autoclave up to 135°C/2 bar (275°F/29 psi)



UNI EN ISO 9001

UNI CEI EN 46001

DENTAL DYNAMIC MANUFACTURING



TeKne Dental s.r.l.

Via del Pescinale, 77 - 50041 Calenzano (FI) - ITALY
+39 055 8825741 +39 055 8825764
info@teknedental.com www.teknedental.com

Distributed by: